



User Manual

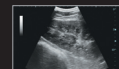
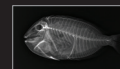




Table of Contents

Important Safety Information	1
Chapter 1: Opening TruDR™	3
Chapter 2: Patient Information	4
Searching and Adding Patients	4
Editing an Existing Patient's Information	6
Creating a New Study for an Existing Patient	7
Chapter 3: Creating a Study Using Modality Worklist.....	8
Chapter 4: Creating Radiographs	11
Choosing Anatomy and Positioning.....	11
Acquiring Radiographs	13
Cropping a Radiograph	17
Chapter 5: Managing Studies	20
Adding New Images to an Existing Study.....	20
Reviewing a Study	22
Burning a CD	23
Exporting a Study as a JPEG.....	25
Sending a Study to a Server	27
Confirming that a Study was Sent	28
Sending a Study to AIS	29
Appendix A: Technical Support.....	30

Important Safety Information

Note



*Read all the instructions before connecting or operating the component
Keep this manual so you can refer to these safety instructions.*

Warning



*There are no user serviceable parts inside.
Refer all servicing to qualified service personnel.*

To reduce the risk of fire or electric shock, do not expose the unit to moisture or water. Do not allow foreign objects to get into the enclosure. If the unit is exposed to moisture, or a foreign object gets into the enclosure, immediately disconnect the power cord from the wall. Take the unit to a qualified service person for inspection and necessary repairs.

Caution



Heed all warnings and safety information in these instructions and on the product itself. Follow all operating instructions.

Clean the enclosure only with a dry cloth or a vacuum cleaner.

Clean the DR plate with a damp cloth. Use no chemicals or alcohol.

You must allow 10 cm or 4 inches of unobstructed clearance around the unit. Do not place the unit on a bed, sofa, rug, or similar surface that could block the ventilation slots. If the component is placed in a bookcase or cabinet, there must be ventilation of the cabinet to allow proper cooling.

Keep the component away from radiators, heat registers, stoves, or any other appliance that produces heat.

Keep the component away from flammable materials.

Place the unit on a fixed, level surface strong enough to support its weight. Do not place it on a moveable cart that could tip over.

Important Safety Information *(continued)*

Connection



The unit must be connected to a power supply only of the type and voltage specified on the rear panel of the unit.

Connect the component to the power outlet only with the supplied power supply cable or an exact equivalent. Do not modify the supplied cable in any way. Do not attempt to defeat grounding and/or polarization provisions. Do not use extension cords.

Do not route the power cord where it will be crushed, pinched, bent at severe angles, exposed to heat, or damaged in any way. Pay particular attention to the power cord at the plug and where it exits the back of the unit.

The power cord should be unplugged from the wall outlet if the unit is to be left unused for a long period of time.

When Problem Occurs



Immediately stop using the component and have it inspected and/or serviced by a qualified service agency if:

- The power supply cord or plug has been damaged.*
- Objects have fallen or liquid has been spilled into the unit.*
- The unit has been exposed to rain.*
- The unit shows signs of improper operation*
- The unit has been dropped or damaged in any way*



CAUTION



**Risk of electric shock.
Do not open.**

Caution: To reduce the risk of electric shock, do not remove cover. No user-serviceable parts inside. Refer servicing to qualified service personnel.



This symbol is to alert the user to the presence of uninsulated dangerous voltages inside the product's enclosure that may constitute a risk of electric shock.

**Applicable for USA, Canada or
where approved for the usage.**


Caution: To prevent electric shock, match wide blade of plug to wide slot. Insert fully.



This symbol is to alert the user to important operating and maintenance (service) instructions in this manual and literature accompanying the product.

Chapter 1: Opening TruDR™

Open the TruDR™ software to begin the process of taking a radiograph.

1. Double-click the TruDR™ icon  on the desktop. This will load your TruDR™ software for acquiring radiographs.
2. Enter your user name and password when the login screen appears.


Your default login is:

User Name: tech

Password: tech



The software will load directly into the Patient Screen if the “Remember Me” box is checked.

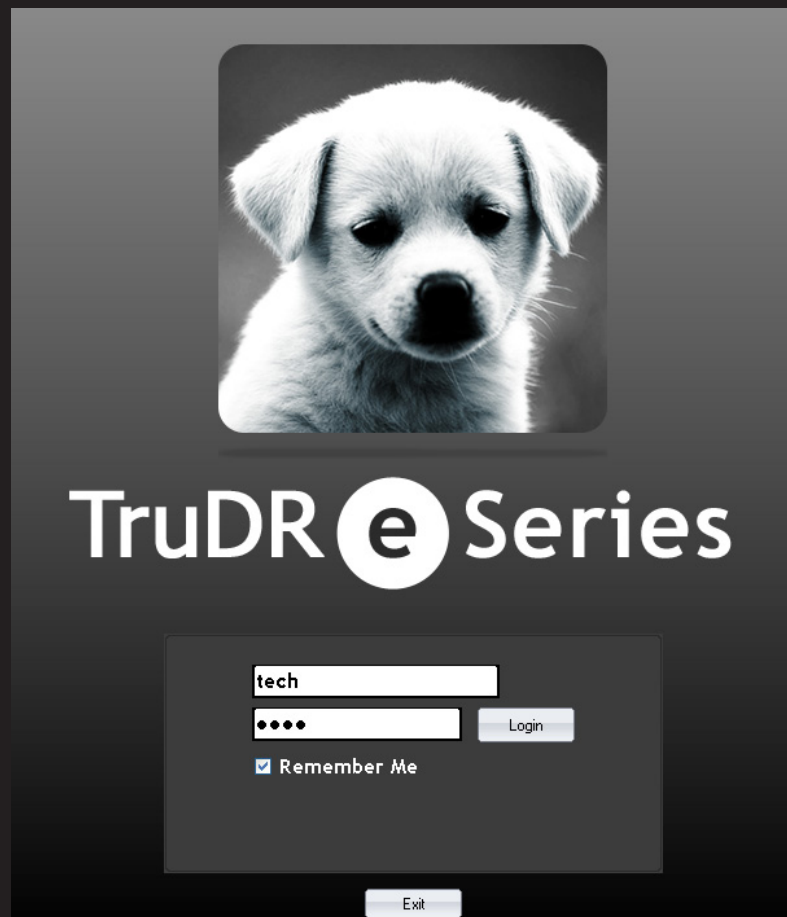

Chapter 2: Patient Information


Searching and Adding Patients


Search for existing patients or add new patients easily by using the filtered search fields.

1. Search for your patient in the Patient Screen by entering the **Patient ID**, the **Patient Name**, or the **Owner Last Name**. Your patient list will automatically filter as you type.



Patients entered within the last 24 hours are displayed in the Patient List.

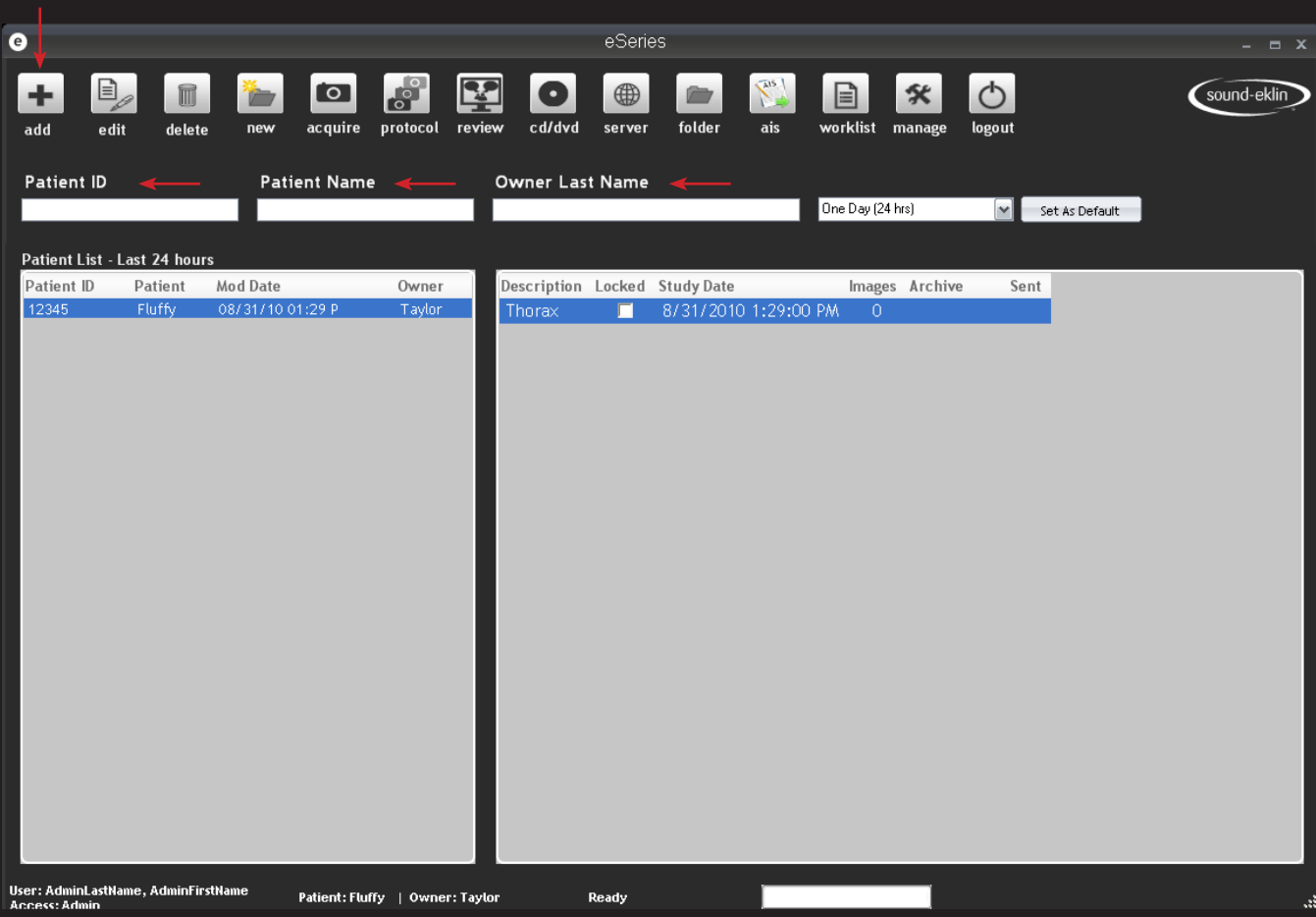


2. If you do not see the patient you are searching for, click the **Add** button  or press the Enter key on the keyboard to bring up the Add Patient Screen.



Chapter 2: Patient Information

Searching and Adding Patients *(continued)*

- 3. Enter your patient's information. **Patient ID**, **Patient Name**, **Species**, **Weight**, and **Last Name** are required fields.
- 4. Click one of the **Save** buttons.

Save		Saves the current patient and returns you to the Patient Screen
Save + Add		Saves the current patient and creates a new patient
Save + Acquire		Saves the patient and prepares the computer to acquire a radiograph.

patient save save+add save+acquire save+protocol

Owners

- James
- Taylor

Patient ID:*

12345678

Patient Name:*

Fido

Secondary Patient ID:

Date Of Birth:

12/03/1998

Species:*

Canine

Breed:

AUSTRALIAN SHEPHERD

Sex:

Spayed

Weight:*

30

lb

Convert to kg

Owner Information

Last Name:*

Dryden

Clear Owner

First Name:

Ron

Institution Name:

Address:

Address2:

Address3:

City:

User: AdminLastName, AdminFirstName
Access: Admin

Patient: Not Selected | Owner: Not Selected


Ready

Chapter 2: Patient Information *(continued)*


Editing an Existing Patient's Information

Edit an existing patient's information for future DICOM tagging.

1. Select your patient to edit from the Patient Screen.
2. Click the **Edit** button. 

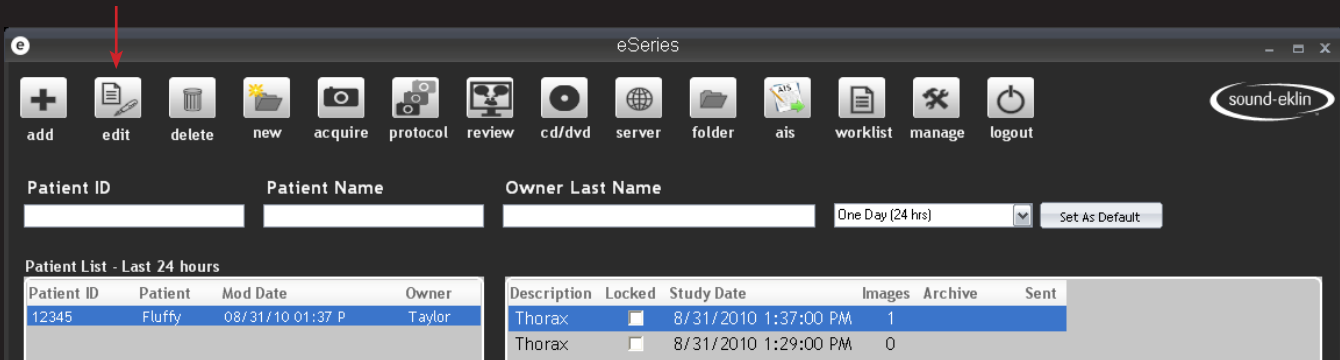

3. Edit the information as needed.
4. Click an appropriate **Save** button (*see page 5*).

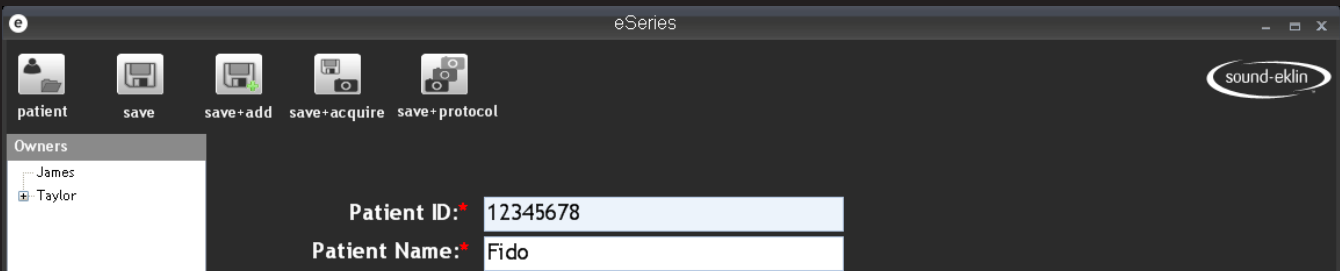



Chapter 2: Patient Information *(continued)*


Creating a New Study for an Existing Patient

Create a new study in a patient previously entered in TruDR™.

1. Search for your patient in the Patient Screen by entering the **Patient ID**, the **Patient Name**, or the **Owner Last Name**.
The Patient List will automatically filter as you type.

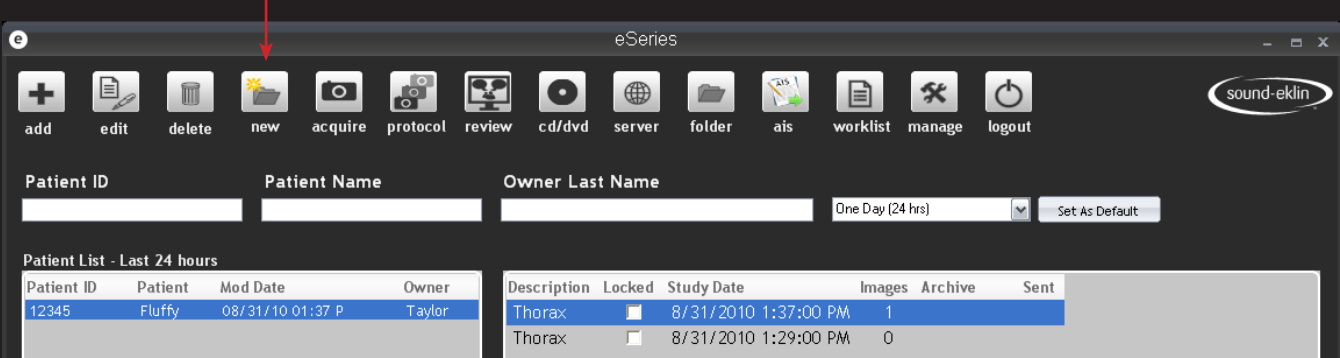


2. Select your patient by clicking on the name in the **Patient List**.
3. Click the **New** button  to create a new study.
4. Select the anatomy and the view from the Anatomy Screen *(see Chapter 4)*.

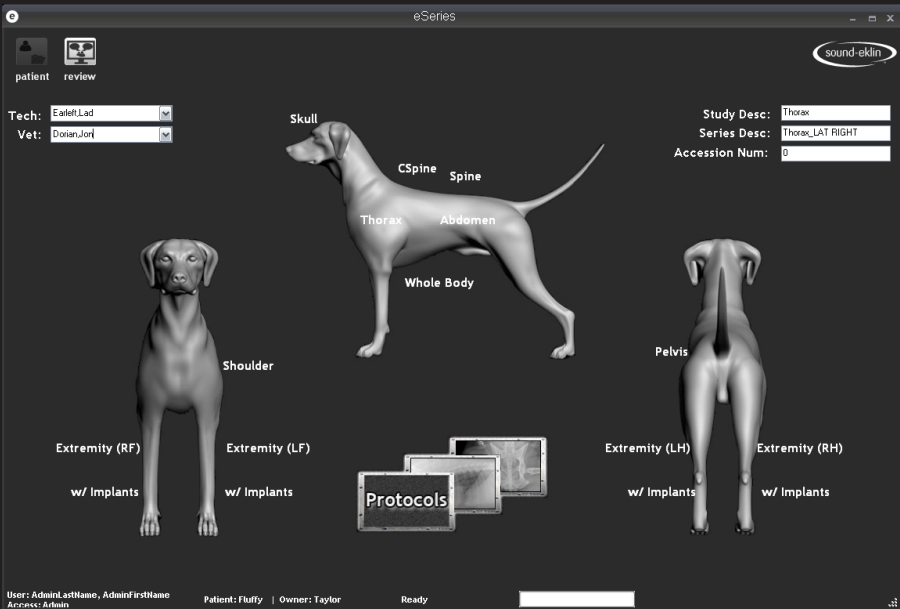


Chapter 3:

Creating a Study Using Modality Worklist

Use Practice Management Software to send the information via Modality.

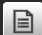1. Use your Practice Management Software to create a modality worklist.
2. In TruDR™, click the **Worklist** button  on the Patient Screen.

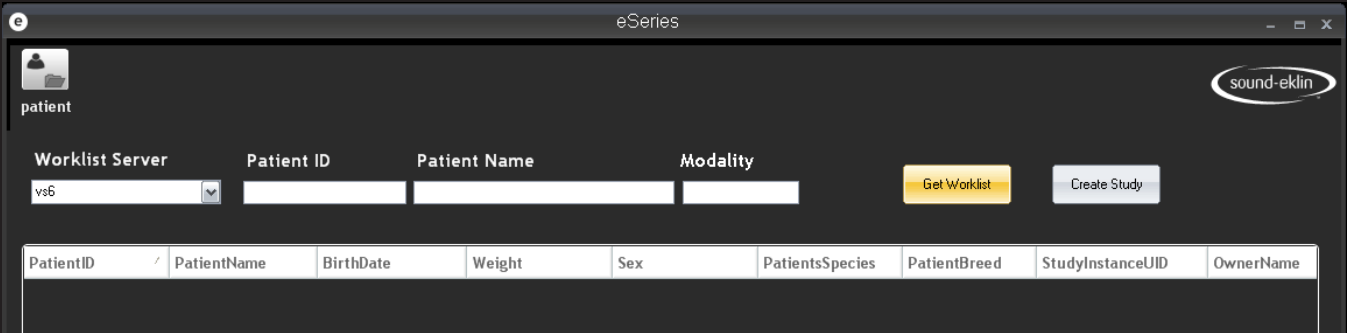


patient


Worklist Server Patient ID Patient Name Modality

vs6

Get Worklist Create Study

PatientID	PatientName	BirthDate	Weight	Sex	PatientsSpecies	PatientBreed	StudyInstanceUID	OwnerName
-----------	-------------	-----------	--------	-----	-----------------	--------------	------------------	-----------

3. Click the **Get Worklist** button to populate a list of all patients sent from your Practice Management Software. You can also search for a patient by entering the patient's information and clicking the **Get Worklist** button. 
4. Select your patient from the populated list.

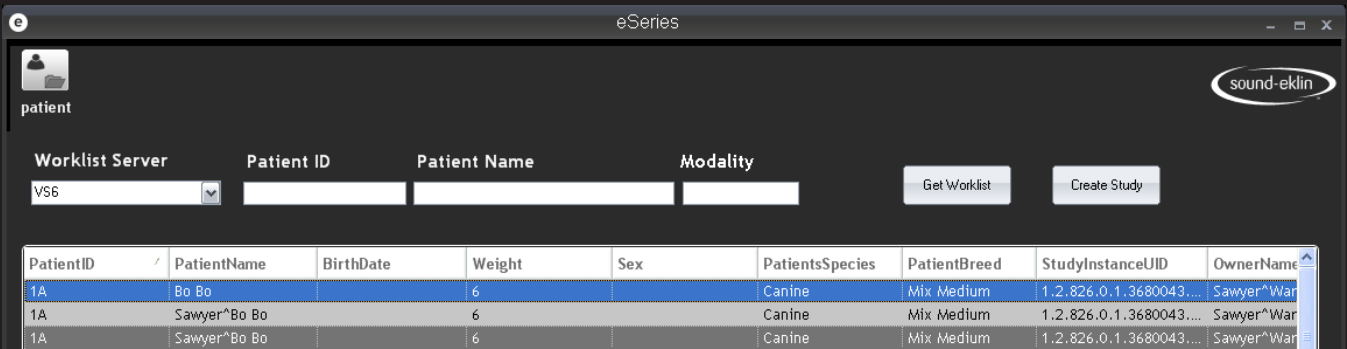

patient

Worklist Server Patient ID Patient Name Modality

VS6

Get Worklist Create Study

PatientID	PatientName	BirthDate	Weight	Sex	PatientsSpecies	PatientBreed	StudyInstanceUID	OwnerName
1A	Bo Bo		6		Canine	Mix Medium	1.2.826.0.1.3680043....	Sawyer^War
1A	Sawyer^Bo Bo		6		Canine	Mix Medium	1.2.826.0.1.3680043....	Sawyer^War
1A	Sawyer^Bo Bo		6		Canine	Mix Medium	1.2.826.0.1.3680043....	Sawyer^War

Chapter 3: Creating a Study Using Modality Worklist *(continued)*


The screenshot shows the 'eSeries' application window. On the left, there's a 'patient' button with a person icon. Below it, there are four input fields: 'Worklist Server' (containing 'VS6'), 'Patient ID' (empty), 'Patient Name' (empty), and 'Modality' (empty). To the right of these fields are two buttons: 'Get Worklist' and 'Create Study'. A red arrow points to the 'Create Study' button. In the top right corner, there is a 'sound-eklin' logo.

- 5. Click the **Create Study** button.
- 6. Fill in any additional information for your patient in the Work List Create Study dialogue box.

The 'WorkListCreateStudy' dialog box is shown. It contains the following fields and values: 'Patient ID' is '1A'; 'Patient Name' is 'Bo Bo'; 'Species' is 'Canine'; 'Breeds' is 'AFFENPINSCHER'; 'Study Description' is 'step desc'; 'Study Types' is 'Abdomen'; 'Tech' is 'TechLastName, TechFir'; 'Vet' is 'VetLastName, VetFirsth'; 'Weight' is '13.2277357310927' lb, with a 'Convert to kg' button next to it. At the bottom, there are 'Save' and 'Cancel' buttons.

- 7. Click the **Save** button.
- 8. Click the **Patient** button  to return to the Patient Screen.

Chapter 3: Creating a Study Using Modality Worklist (continued)

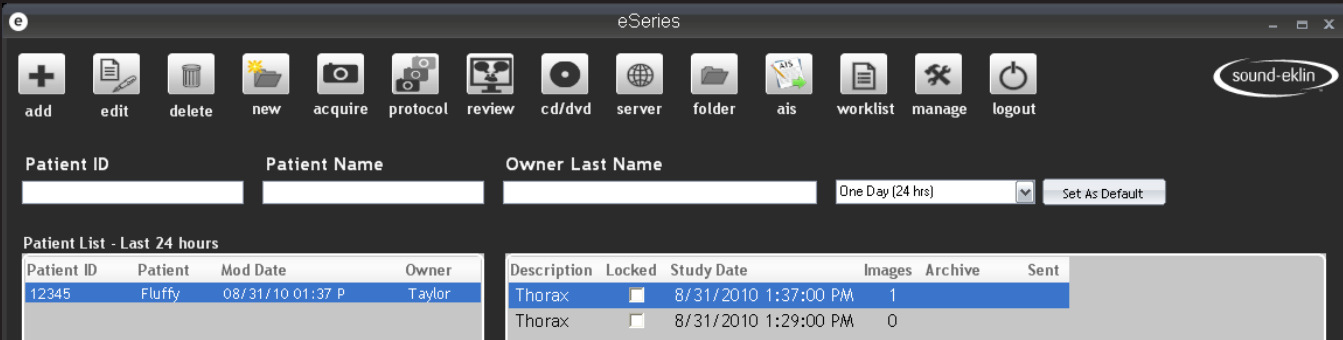


9. Select your patient and their study from the Patient Screen.



The number of images in the study reads “0” since no images have been captured.

10. Click the Acquire button to begin acquiring radiographs (see pages 12-14).

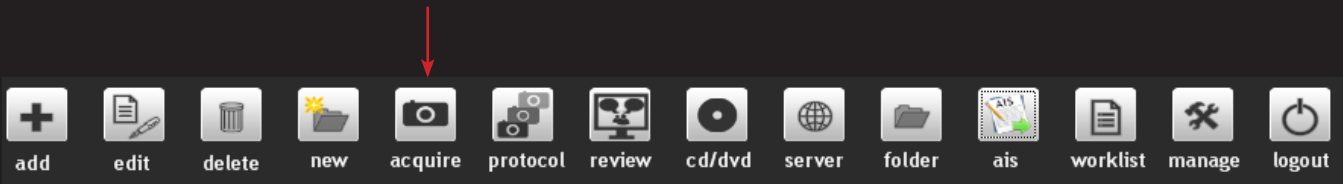


Chapter 4: Creating Radiographs

Choosing Anatomy and Positioning

Choose the anatomy and positioning to automatically enhance and position your radiograph.

Tech: Earleft,Lad
Vet: Dorian,Jon

1. Select the Tech and ordering Vet from the drop-down menus.
2. Select the appropriate anatomy by clicking on the corresponding text. This will bring you to the Views Screen.

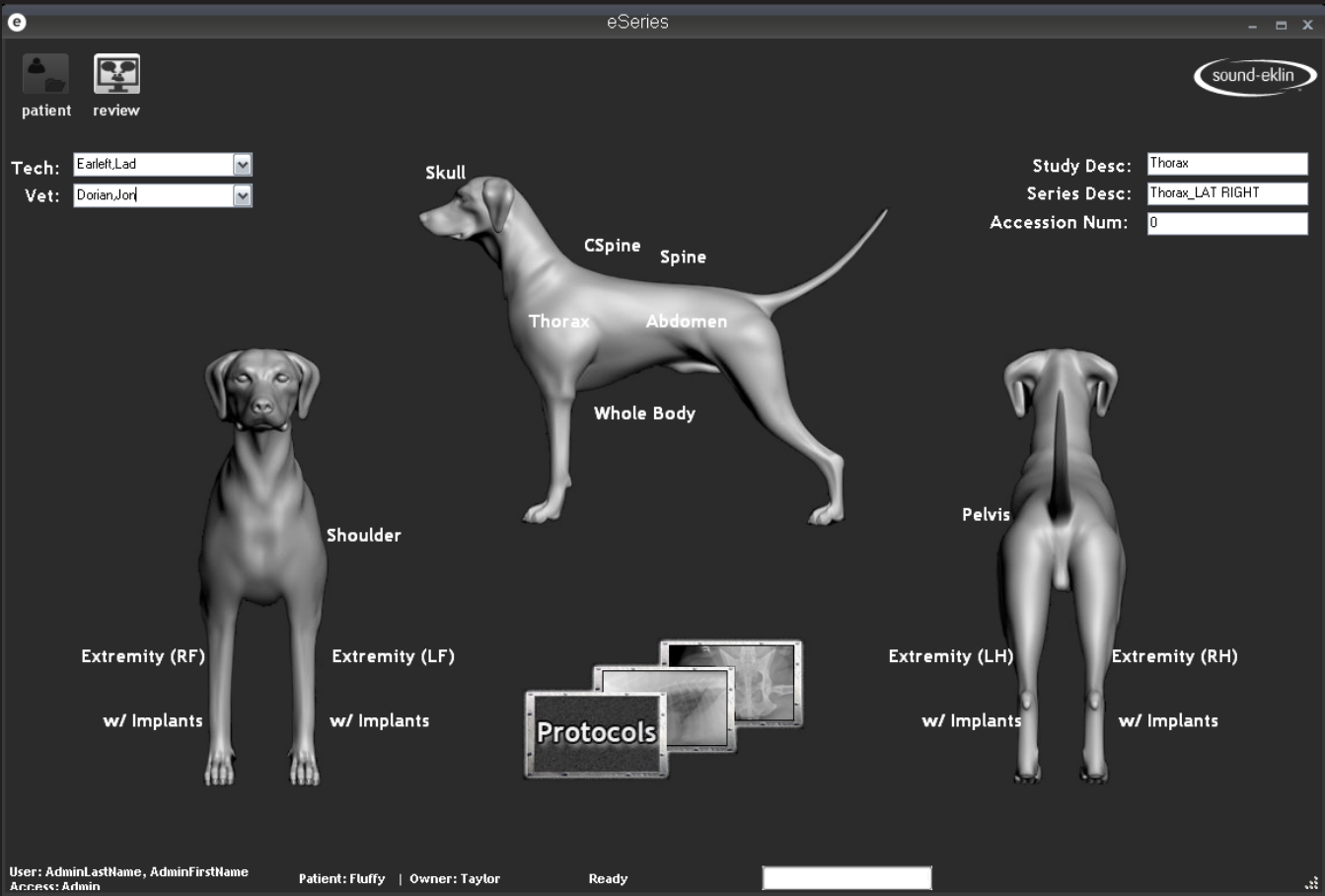

Chapter 4: Creating Radiographs

Choosing Anatomy and Positioning (continued)

3. Click on the appropriate view button for your shot. This will bring you to the acquisition screen.



Chapter 4: Creating Radiographs *(continued)*

Acquiring Radiographs

1. Set the X-ray machine according to the Sound-Eklin™ Technique Chart using the appropriate KVP and mAs.




Due to variations between x-ray generators, optimal techniques may vary somewhat from the values printed on the Sound-Eklin™ Technique Chart.

2. Place the patient on the table with the desired anatomy in the center of the light field cross hairs.
3. Press the foot pedal to prep and to expose the X-ray machine.



Proper Two-Stage Foot Pedal Technique


1. To Prep, press the foot pedal halfway down.
2. Wait until the Generator is ready to fire.
3. To expose, press the foot pedal completely down.

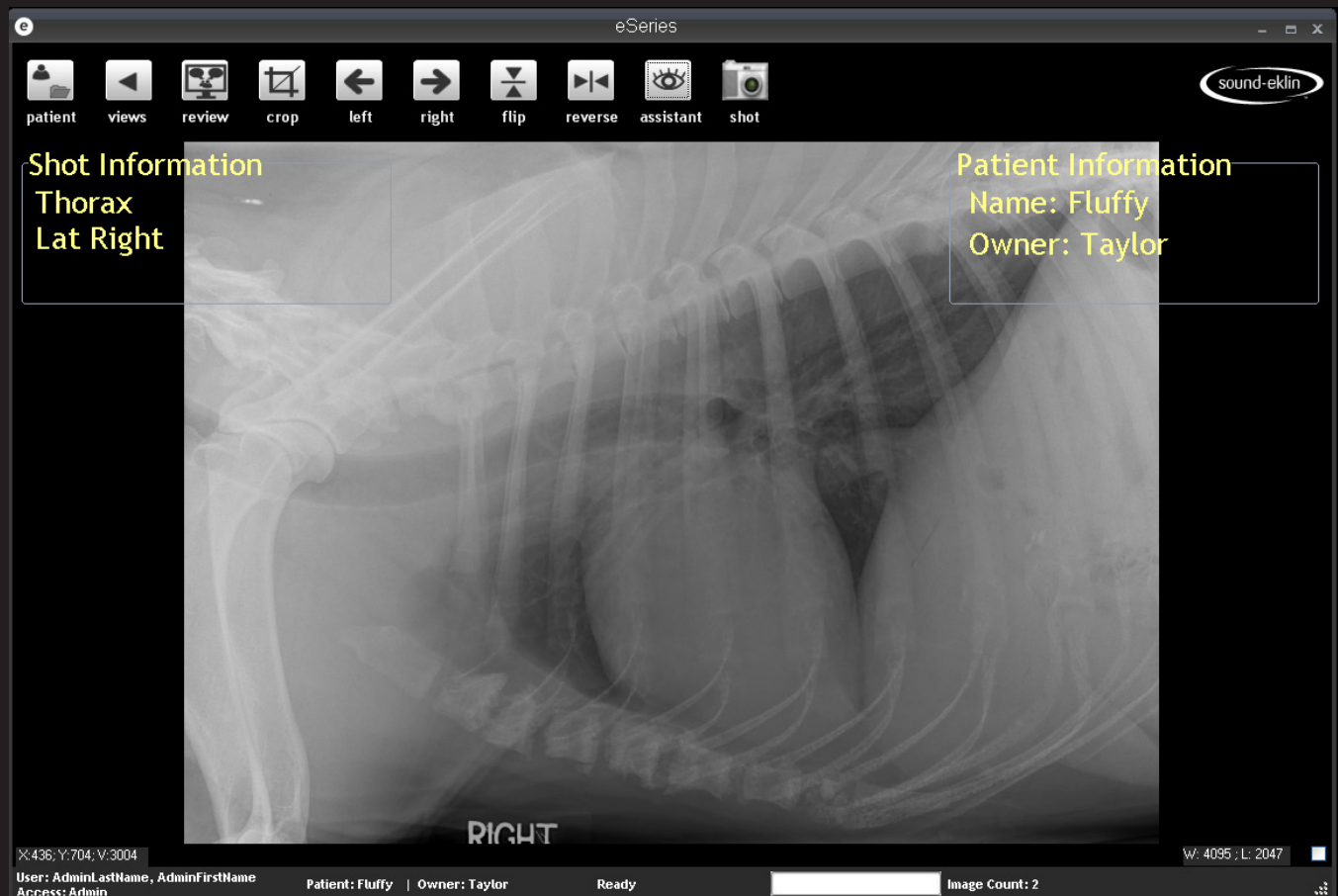


Chapter 4: Creating Radiographs

Acquiring Radiographs *(continued)*

- 4. The image will appear shortly. Use the **Left**, **Right**, **Flip**, and **Reverse** buttons to make orientation adjustments before shooting your next shot.

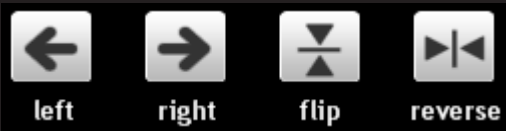

- 5. After image acquisition, the Positioning Assistant displays in the lower right corner of the screen.

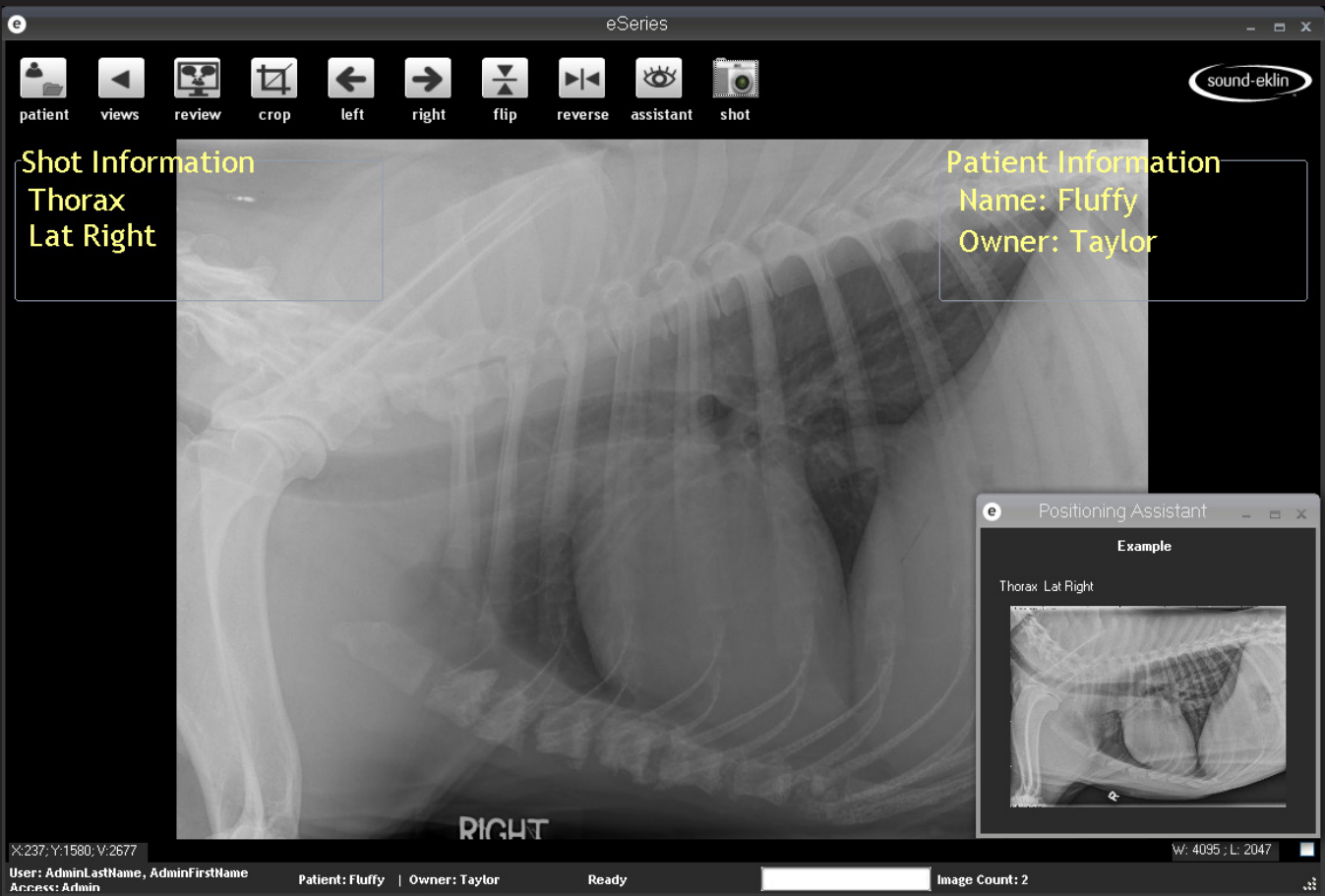

Chapter 4: Creating Radiographs

Acquiring Radiographs *(continued)*



- The pop-up window lists the anatomy and view, and shows the correct position of the animal.



Click and drag to relocate the Positioning Assistant pop up window.

Click the minus icon to minimize the Positioning Assistant, the expand to maximize icon, and the “x” to close.



Select the assistant button to toggle Positioning Assistant off/on (Enable/Disable).





Crop your image as needed (see pages 17-18).

If you need to reshoot your X-ray, press the foot pedal and fire again.

Chapter 4: Creating Radiographs

Acquiring Radiographs *(continued)*

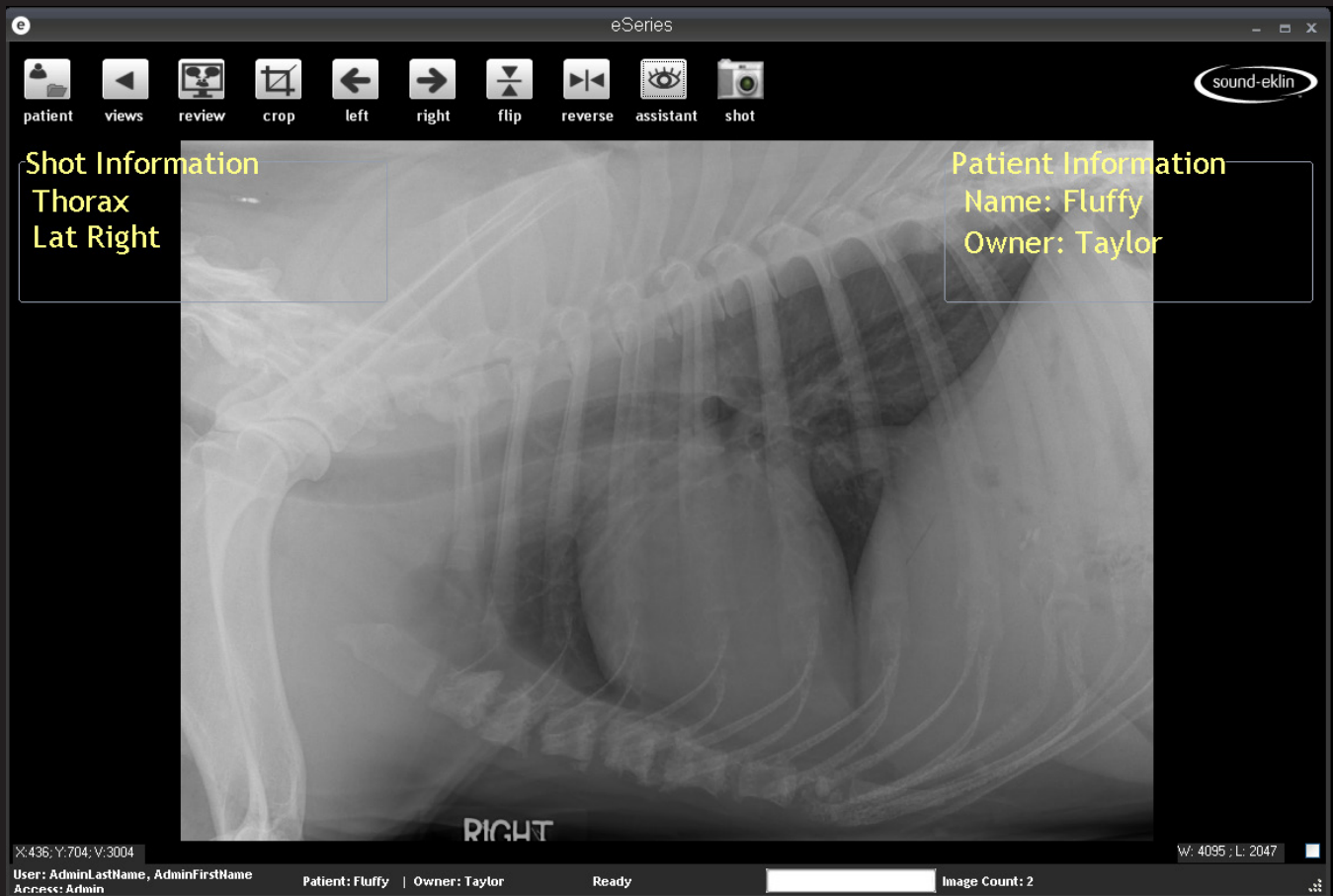


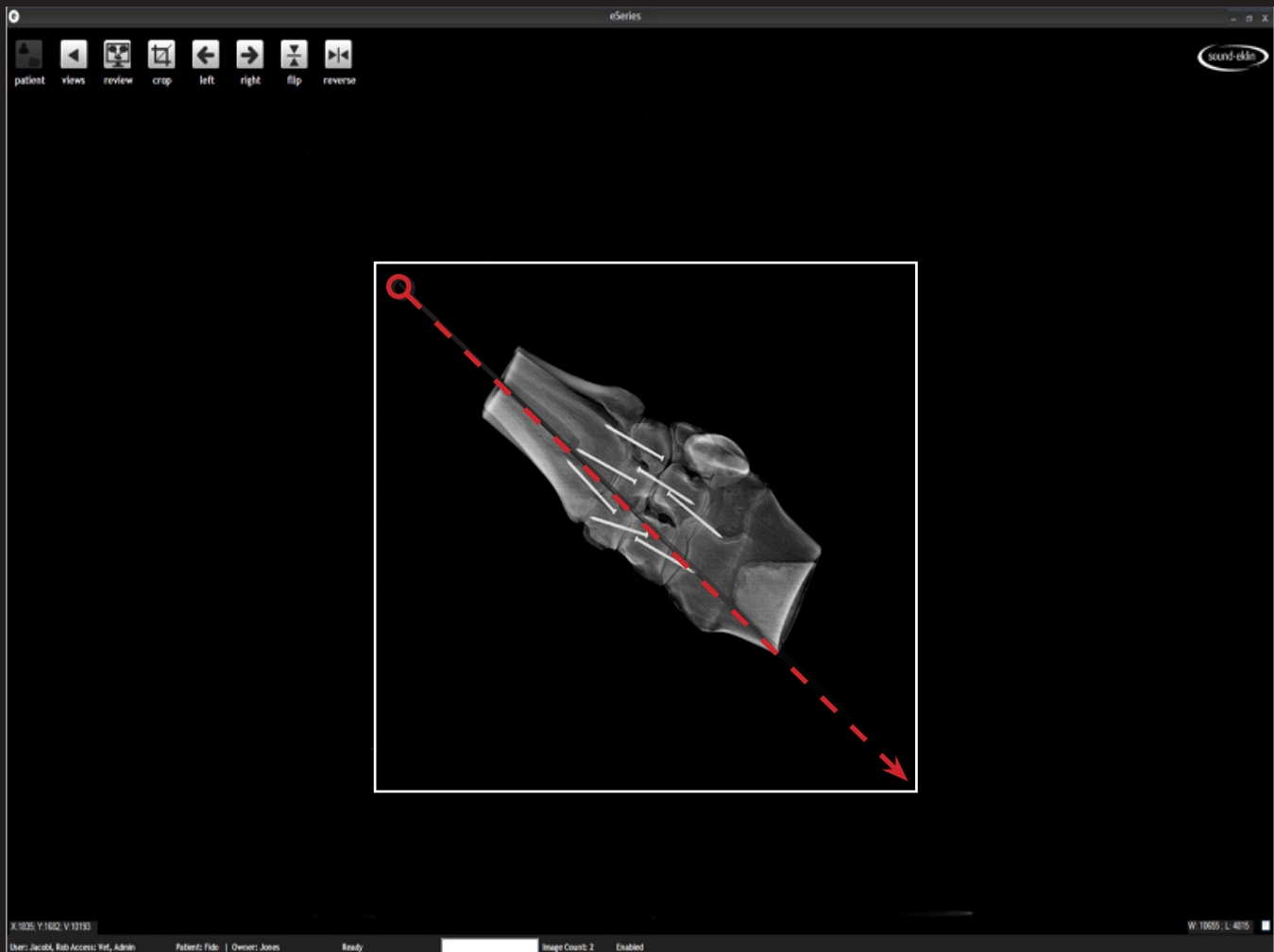
7. Click the **Views** button to change views. 
8. Select your next view from the Views Screen.
Click the **Anatomy** button from the Views Screen if another anatomy is needed *(see page 11 and 12)*.

Chapter 4: Creating Radiographs *(continued)*

Cropping a Radiograph

Crop out unwanted portions of your radiograph.

1. Click the Crop button in the Acquisition Screen. 




The image will enhance once the cropping is complete.

Chapter 4: Creating Radiographs

Cropping a Radiograph (*continued*)


2. Frame the desired image by holding the left mouse button in the upper left corner of the image and dragging to the lower right hand corner.

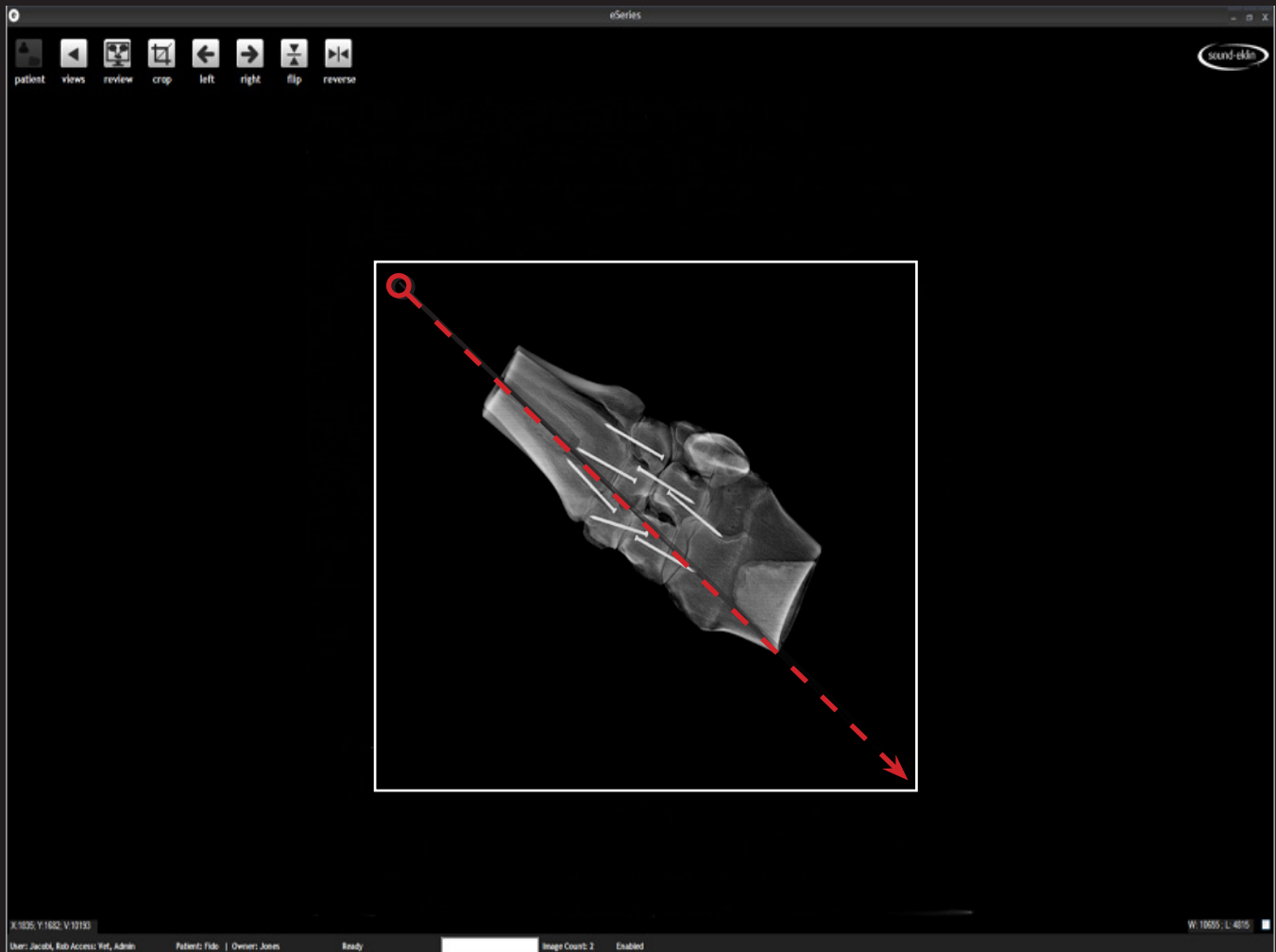

Click and drag your mouse by starting at the circle and dragging down toward the arrow. The arrow will not appear in the actual software.

Chapter 4: Creating Radiographs

Cropping a Radiograph (*continued*)

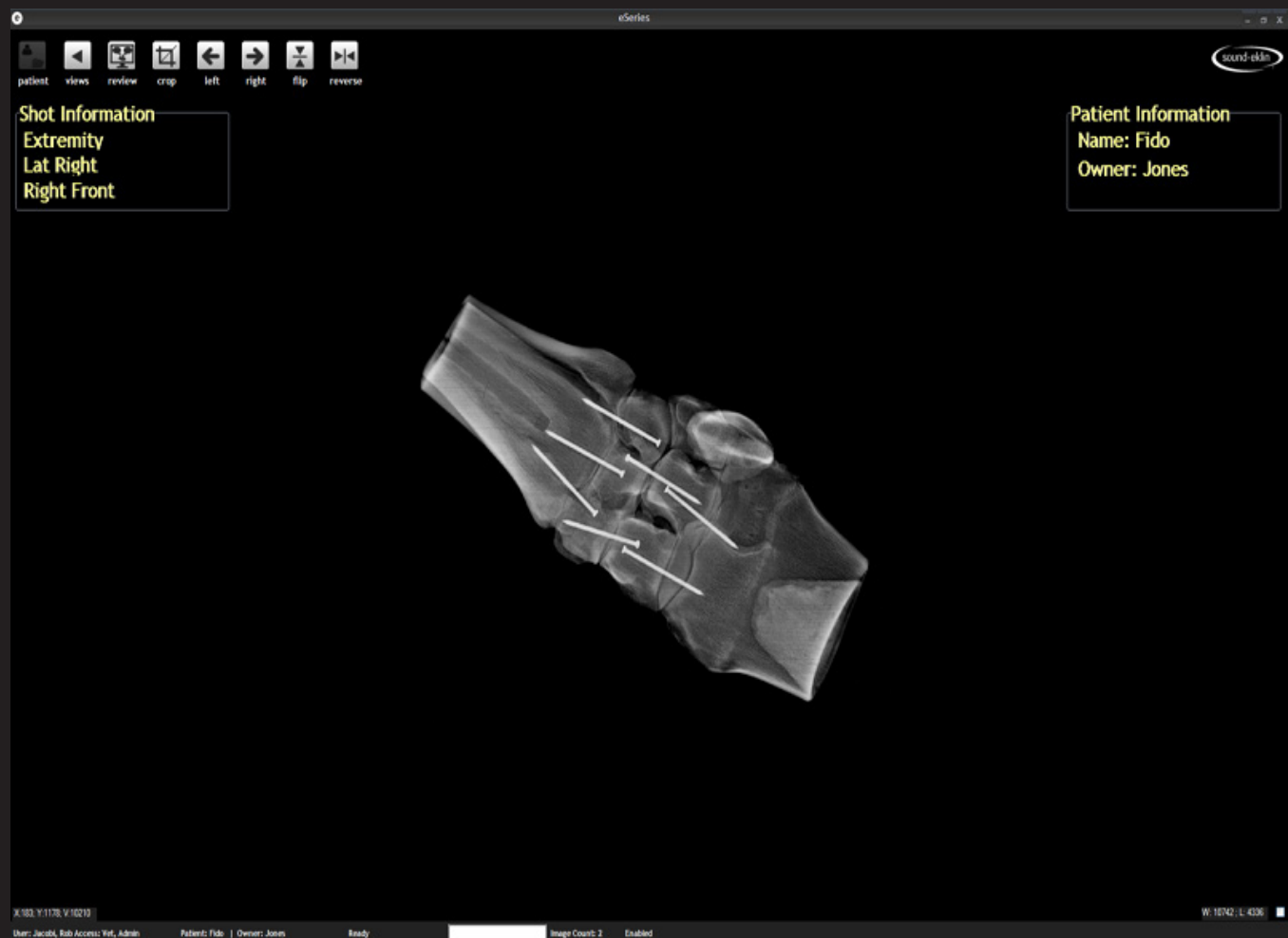

3. The image will enhance once it is cropped.
4. Re-crop your image if necessary before taking your next radiograph.




Images must be cropped before shooting the next image.

Chapter 5: Managing Studies

Adding New Images to an Existing Study (continued)

3. Select a view and acquire your images (see page 13).

 Previous shots will load in the thumbnail pane on the left side of the Review Screen. New shots will appear at the bottom.






Chapter 5: Managing Studies *(continued)*

Reviewing a Study


Change the layout or delete any unwanted images before your study is archived on the Server.

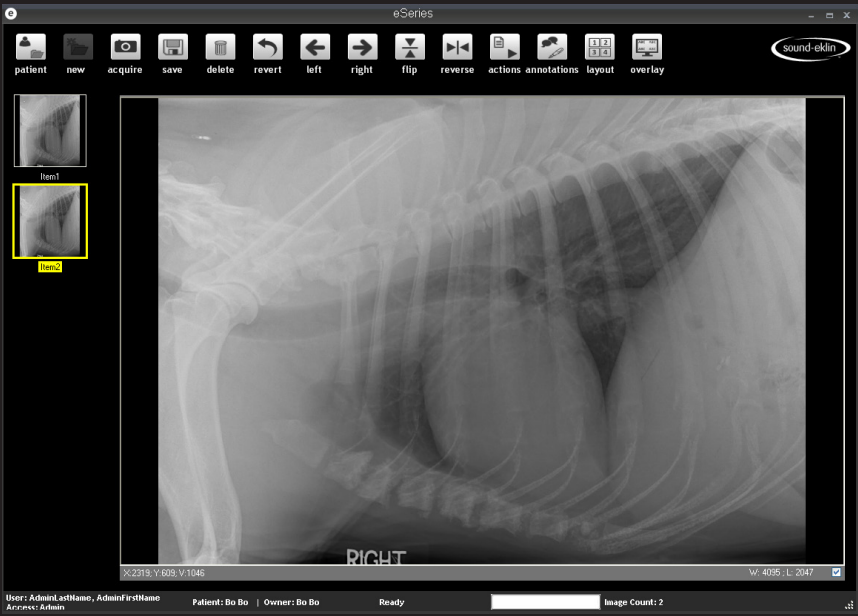1. Select a thumbnail in the left-hand pane to load the image into the review pane on the right. The selected image is highlighted by a yellow box.



2. Click the Delete button to remove any unwanted images. 
  *Deleting an image permanently removes it from the computer.*

3. Review the images by using the following techniques:
 Zoom: Roll the mouse wheel forward and back.
 Pan: Click and hold the left mouse button and move the mouse.
 Window Level: Click and hold the right mouse button, then move the mouse diagonally across the image.


4. Click the Patient button  to end your study and transfer the X-rays to a Review Station/Server. The software will return you to the Patient Screen.


Chapter 5: Managing Studies (continued)

Burning a CD

Burn a copy of a study to a CD.

1. Search for your patient in the Patient Screen by entering the **Patient ID**, the **Patient Name**, or the **Owner Last Name**. The Patient List will automatically filter as you type.
2. Select your patient and your study. Your patient and study are highlighted in a blue box.

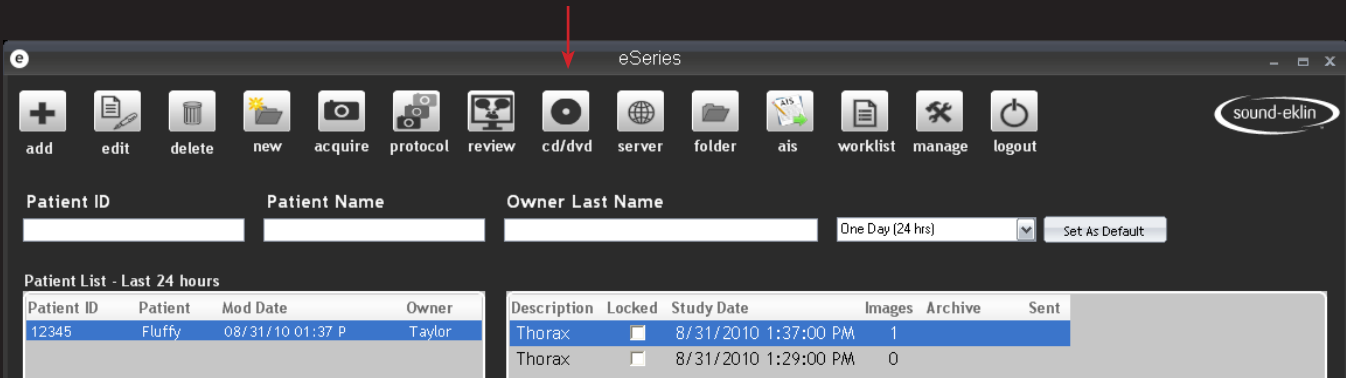

3. Click the **CD/DVD** button. 
4. Select an image format from the **Export As** drop-down menu.



DICOMDIR
JPEG

DICOM images require a DICOM viewer
JPEG images can be viewed on any computer

5. Click the **Add DICOM Viewer** check box if you selected the **DICOM/DIR** option from Step 4 above.



The options are unchecked by default

Annotation Burning

Imprints annotations made to saved images

Anonymize Images

Strips all client and clinic information from the images

Repository Format

Places images in an equestrian repository folder format

Add DICOM Viewer

Adds a DICOM viewer to view the images on a computer

6. Click the **Select** button.


Chapter 5: Managing Studies *(continued)*

Burning a CD *(continued)*

7. Insert a blank CD (CD-R or CD-RW) into the computer.
8. Click the **Write Disc** button.



*The CD burn progress is shown in the lower right hand corner.
You will be prompted to remove the disc when the burn is finished.*

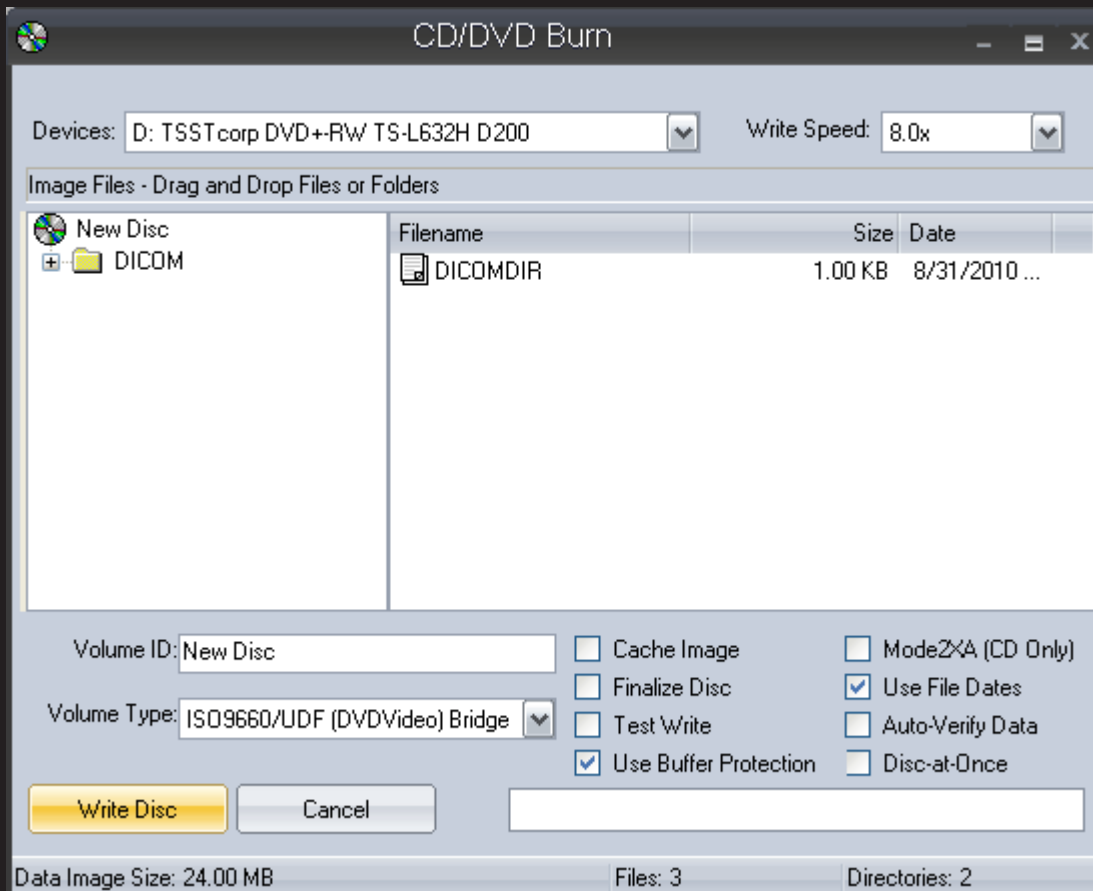


9. Click the **OK** button to complete writing the CD and remove the CD from the CD-Rom.

Chapter 5: Managing Studies (continued)

Exporting a Study as a JPEG


Export images from a study to a JPEG format.

- 1. Select your patient and your study to export from the Patient Screen.
- 2. Click the **Folder** button. 
- 3. Select the JPEG image format from the **Export As** drop-down menu.
- 4. Select an **Overlay Burning** option from the drop-down list.



The default option is set to None

- None** Shows no overlays on JPEG image
- Detail** Shows a detailed DICOM overlay on the JPEG image (shows image number, anatomy, technique, hospital, vet, patient, client, species, DOB, sex, acquisition date)
- Summary** Shows a summary DICOM overlay on the JPEG image (shows anatomy, view, laterality, hospital, vet, patient, client, species, DOB, sex, acquisition date)

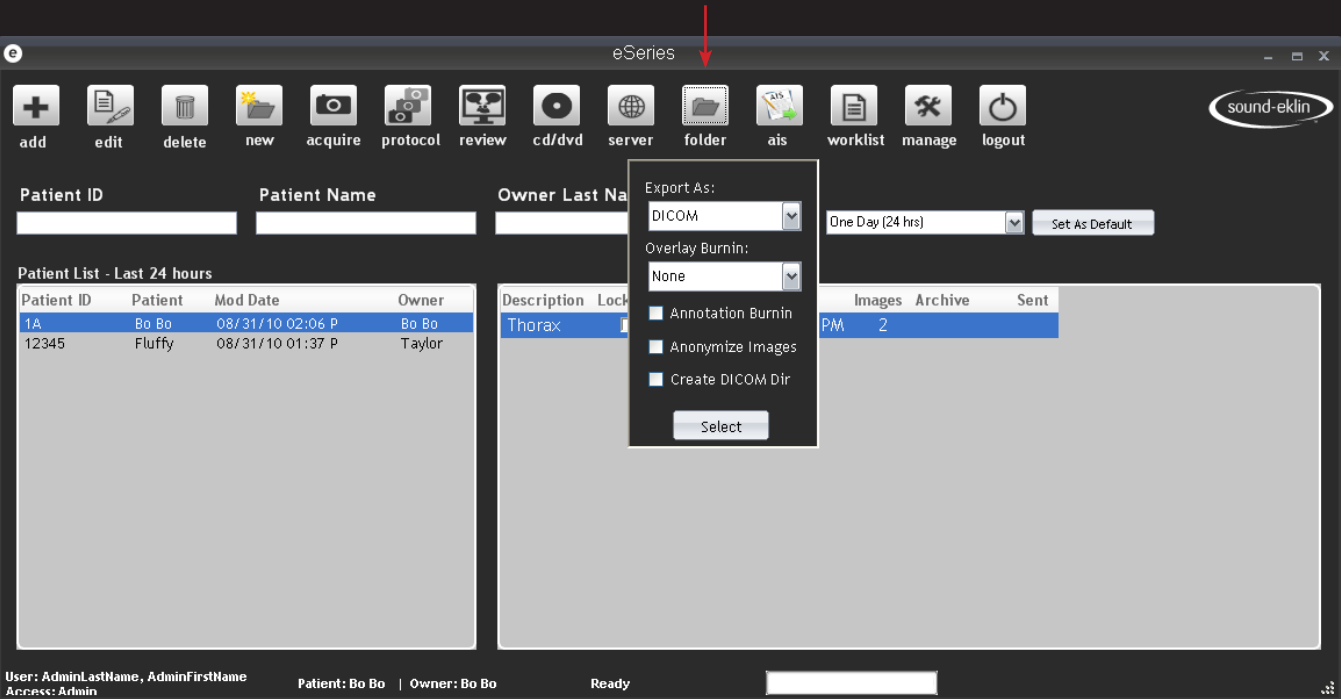

Chapter 5: Managing Studies

Exporting a Study as a JPEG (continued)

5. Check any additional options boxes you want to add to your images.



The options are unchecked by default

Annotation Burning *Imprints annotations made to saved images*

Anonymize Images *Strips all client and clinic information from the images*


Repository Format *Places images in an equestrian repository folder format*

6. Click the Select button.

7. Select a location to save your images in.




Your images are saved to a selected location in a folder with the animals name and the clients last name (e.g. Fido-Smith).

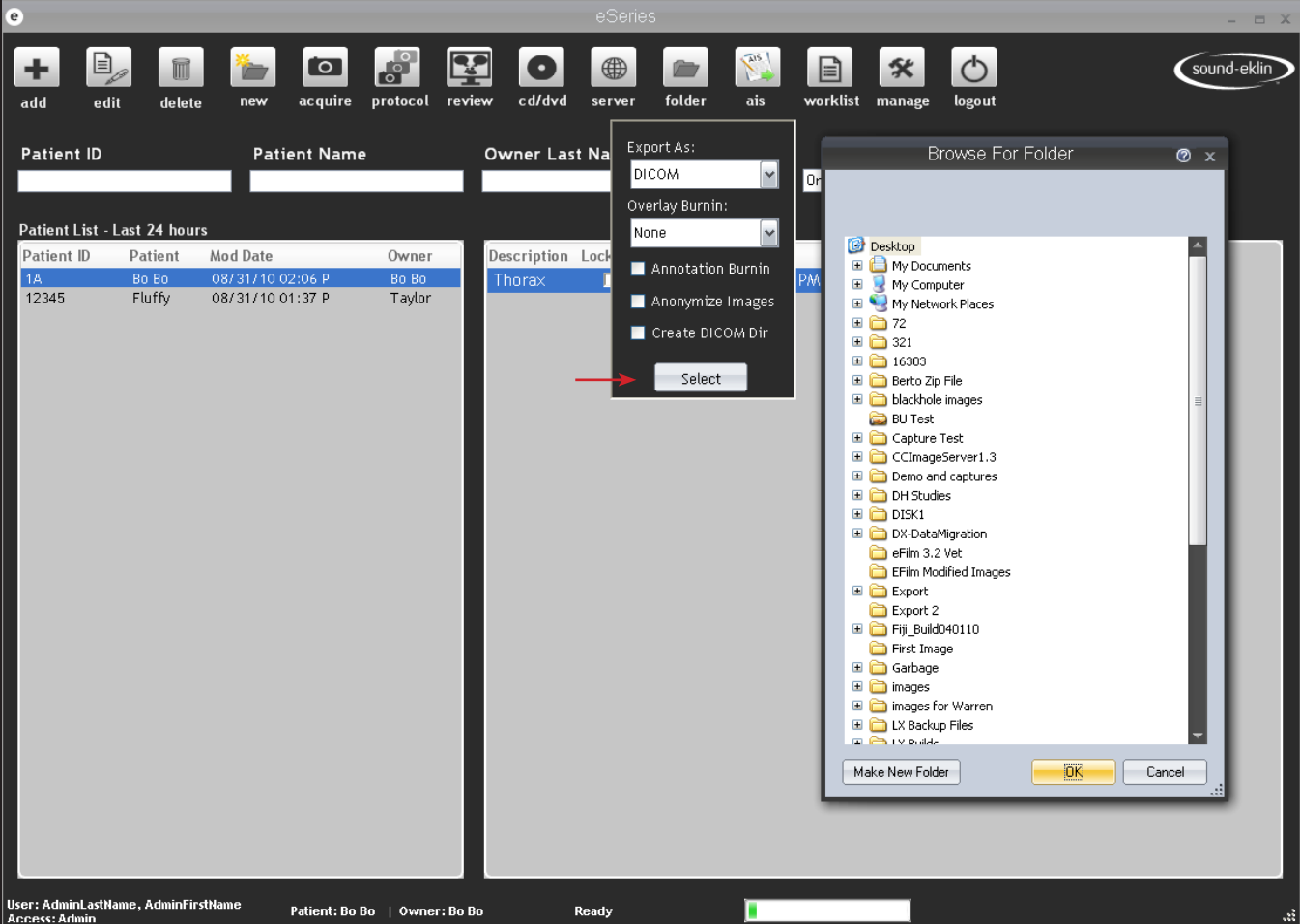



Chapter 5: Managing Studies (continued)

Sending a Study to a Server

Send a study to a DICOM server.

- 1. Search for your patient's study by entering the **Patient ID**, the **Patient Name**, or the **Owner Last Name**. The text will automatically filter as you type.
- 2. Select your patient by clicking on the name in the **Patient List** field.
- 3. Click the **Server** button. 

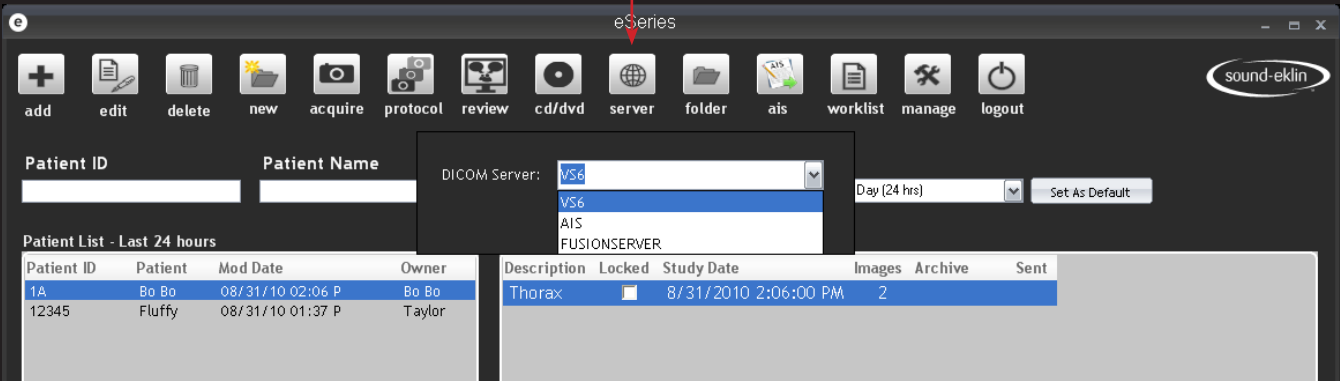

- 4. Select a **DICOM Server** from the drop-down list.
-  *DICOM server information will be provided by the recipient.*
- 5. Click **Select**.
- 6. Click the **Start Store** button to begin the DICOM image transfer.




Confirming that a Study was Sent

1. Finish your study by returning to the Patient Screen.
The Archive field will show a server name (e.g. SOUNDSERV1).
2. Hold the mouse over the server name to show a pop-up window with the transfer data.




Green Box	<i>Image transfer complete</i>
Red Box	<i>Image transfer incomplete or failed</i>

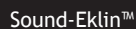

Chapter 5: Managing Studies (continued)

Sending a Study to AIS

Send a telemedicine study to AIS.

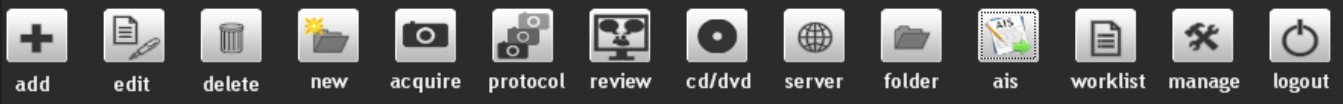

1. Select the patient you want from the **Patient List** field and click the **AIS** button from the toolbar.



Once you have clicked the AIS button the images are automatically sent to AIS and a New Consultation window is launched.

2. Fill in all the required information on the **New Consultation** window and click **Submit**.

New Consultation

Consultation Information

Selected Practitioner * : Vetfirstname Vlastname of Sound

Patient Name *

cat

Owner Last Name *

Henderson

Owner First Name *

Warren

Species *

Feline

Breed

Bengal

Other Breed (if not in breed list)

Date of Birth *

1999

September

19

Color

Sex *

Female

Chip No.

Case No.

cat

Owner Phone

Priority & Response Time

24

hrs

Select Specialty *

Select specialties by moving a highlighted specialty from left to right. Order the specialties by using the up-down carats below the right select field.

Available Specialties

AAHA Certification

Cardiology

Cardiology with Verbal

Clinical Pathology

Computed Tomography

Dentistry

Dermatology

ECG - Event Monitor

ECG - Holter Monitor

ECG - Preop

ECG - Routine

ECG - Stat

Magnetic Resonance Imaging

Medicine

Medicine with Verbal

Current Specialties

Radiology

Wt: 0 lbs

Temperature: C

Pulse: bpm

Respiration: per minute

History (Not required for Petnet®)

Description of the Case

Practitioner Comments (private comments to the specialist - these will NOT appear in the report)

Submit

Cancel Consultation

Appendix A: Technical Support

Shipping Address

Sound-Eklin™

5817 Dryden Place, Suite 101

Carlsbad, California 92008

USA

Technical Support

Office hours: Weekdays 8:00 A.M. — 5:00 P.M. Pacific Time

Emergency 24-hour support is available.

Telephone: 800-268-5354, option 3

Website

www.soundeklin.com



Sound-Eklin™ Technical Support

5817 Dryden Place, Suite 101
Carlsbad, California 92008 USA

toll free: 800.268.5354, option 3
telephone: 760.918.9626
fax: 760.918.9620
international: +1.760.918.9626
www.soundeklin.com

© 2010 All Rights Reserved
Sound-Eklin™ — a VCA ANTECH company